



Infusion Solution



Infusion workstation
Central infusion monitoring system

Smart Infusion, Lead Future


Upgrade Large Touch Screen, Easy to Operate Clearly and Lvely

- High-definition screen: 4.2inch IPS HD touch screen, viewing suppor 170° for all data showing
- Sensitive operation: Millisecond level response without any delay for operation even with gloves
- Automatic pump door: open the door with one button in place to reduce the misoperation of door open
- Manual and Auto-integrated push handle: intelligent sense of syringe push handle position, automatically eliminate sensor clearance, fast start infusion



Medical Advice Linked with Drug Library, Double Assurance Infusion Safety


- Medical advice issue: Easy to connect with HIS, automatic issue the medical advice to the pumps to realize the close system

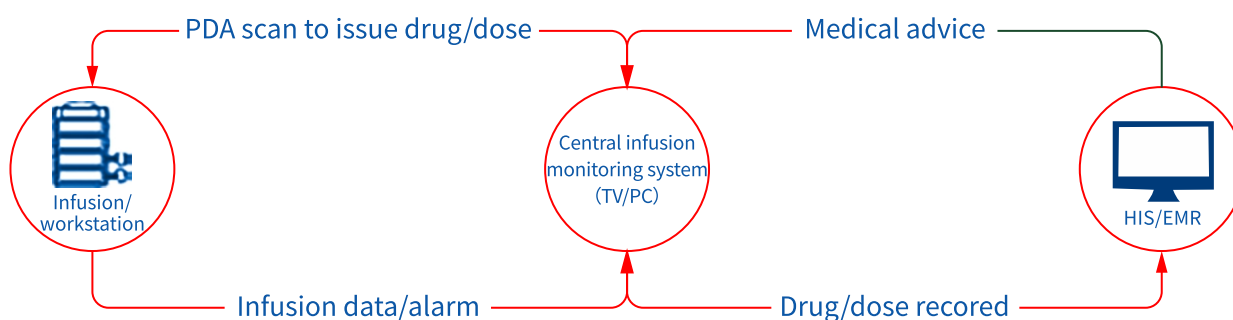

- Drug library support: Equipped with drug color recognition system, built-in drug library and DERS, instant warning of specification and dosage errors




Drug selection			
■	Propofol	10mg/mL	500mg, 50mL
■	Fentanyl	0.02mg/mL	1mg, 50mL
■	ADREnaline	1mg/mL	50mg, 50mL
■	Sufentanil	0.4mg/mL	20mg, 50mL

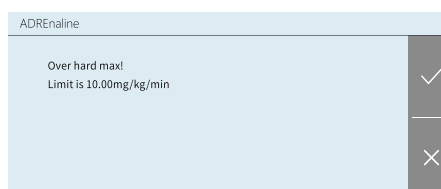

System Interconnected, Information Integrated

- Remote monitoring: Use the central monitor system or mobile terminal to collect the infusion data, achieve the remote monitoring out of the ward
- Nursing record: Automatically record and load the infusion time, dose volume, infused volume to the system, save the time on paper working
- Data sharing: Support HL7 protocol, infusion log data is open to share, collect the data in time, and rejection of information islands

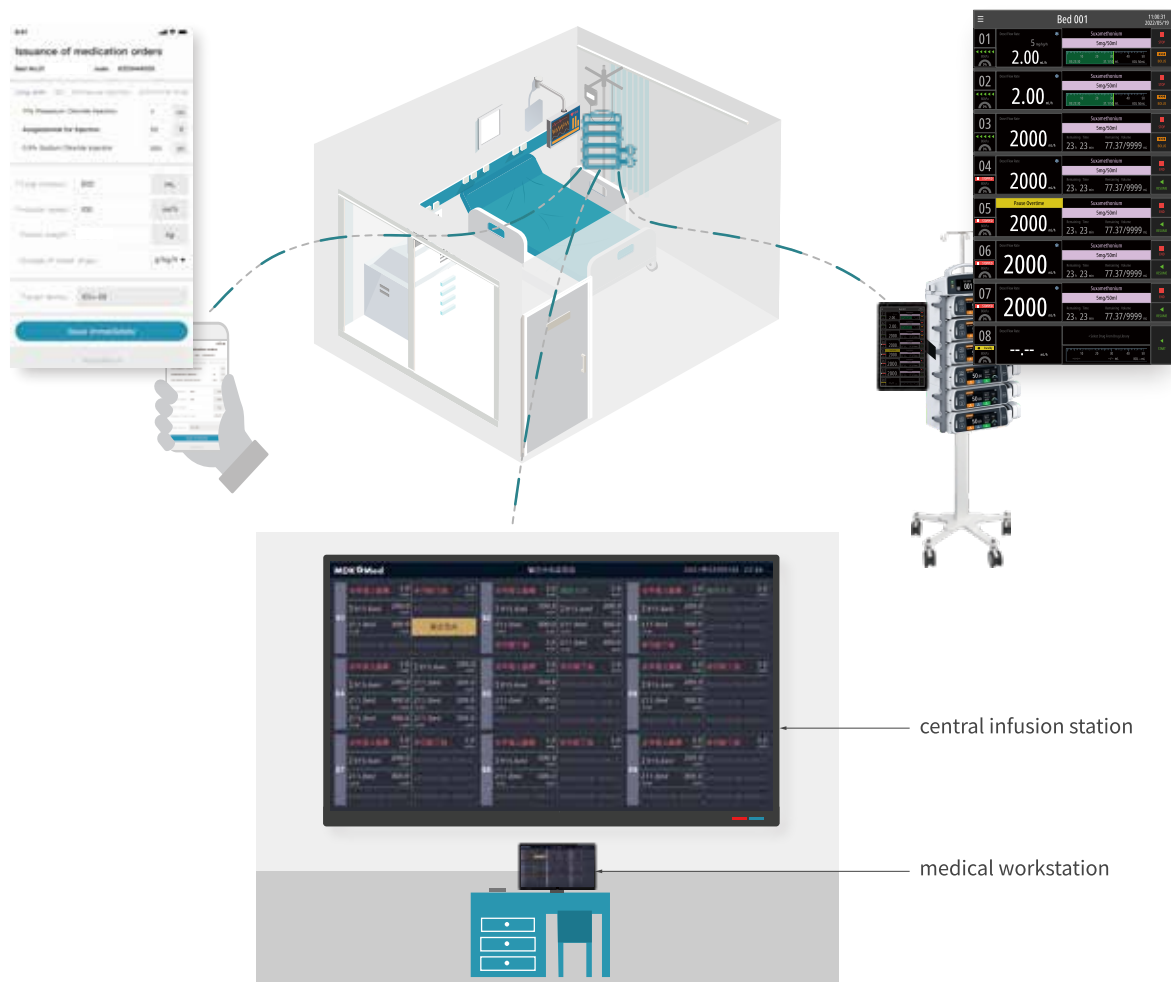

Flexible Docking, Easy to Stack

- Flexible docking: Free combination to 2-24 channels, match various infusion scenes easily
- Easy insertion and removal: the sequence of single pump can be adjusted freely, which is more convenient for clinical application
- Safety warming: Support constant heating during the infusion

2 channels(general ward)



4-8 channels(ICU)



10-24 channels(ICU)



Product parameters

HK series infusion pump



HK-100I

Hardware configuration	4.2-inch HD screen	Automatic pump door
Mode	RVT, Dose, Weight, Relay, Gradient, Drip, First Delivery, Drug, Loading Dose, Micro, Sequence, Intermittent, TCL, PCA, TPN, Blood Transfusion, Oncology Therapy;	
Volume	0.01~9999.99mL, step by 0.01mL	
Rate	0.01-2000mL/h, step by 0.01mL/h	
Sound	≤40dB	
Accuracy	±4%	±3% Dynamic pressure sensitive detection (dedicated consumables)
DERS	DERS intelligent medication system	
Information Transmission	Support WIFI, HL7, Wired	
Dimensions	216×122×72mm(width×depth×height)	1.3kg

HK series syringe pump



HK-100II

Hardware configuration	4.2-inch HD screen	Manual and Auto-integrated push handle
Mode	RVT, Dose, Weight, Relay, Gradient, First Delivery, Drug, Loading Dose, Micro, Sequence, Intermittent, TPN, Blood Transfusion, Oncology Therapy;	
Volume	0.01~9999.99mL, step by 0.01mL	
Rate	0.01-2000mL/h, step by 0.01mL/h	
Sound	≤40dB	
Accuracy	±1.5%	High sensitive detection of occlusion pressure(dedicated extension tube)
On-site detection	Yes	
DERS	DERS intelligent medication system	
Information Transmission	Support WIFI, HL7, Wired	
Dimensions	275 x125×72mm(width×depth×height)	1.5kg



PCA unit



Handle



Data interconnection



Quick mount



Adapter

Medi Waves INC

Office: B-68, G.T. Karnal Road Industrial Area, Delhi-110033, India
Works: 92, Sec-57, Phase-IV, HSIIDC Kundli, Sonapat, Haryana - 131028

✉ info@medi-waves.com
☎ +91-11-4238-4227