




MEDI WAVES INC.

for a healthier tomorrow...



P-003DX

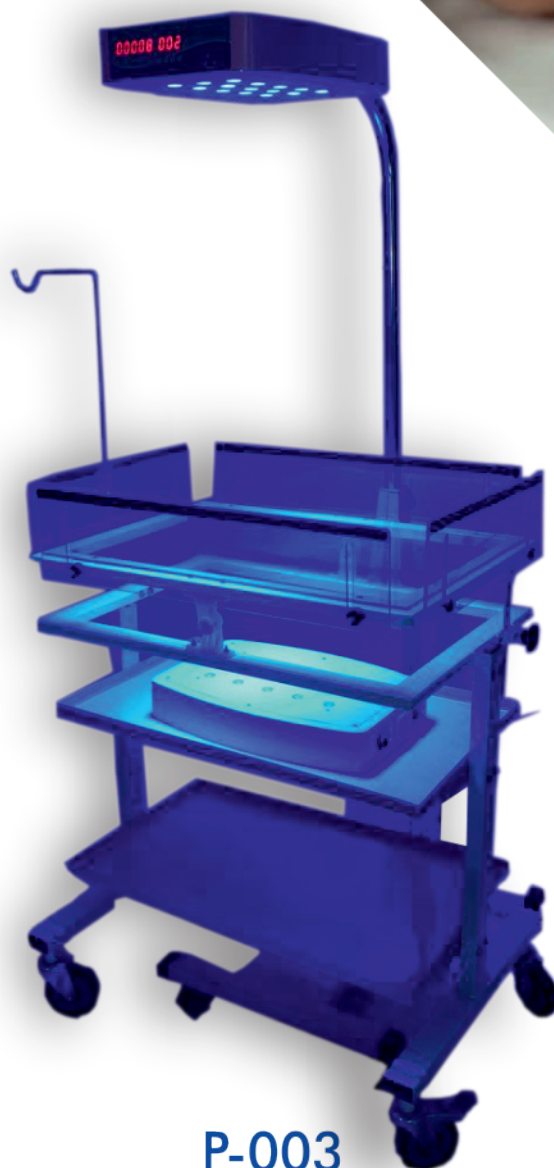


P-003



PHOTOTHERAPY UNIT

LED PHOTO THERAPY UNIT

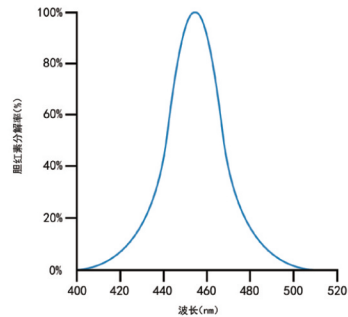

- Phototherapy based on LED technology, which after filtering provide, a light of wavelength approximately 450 to 470 nm with peak wavelength of 450-460nm range.
- Irradiance minimum 35 $\mu\text{W}/\text{cm}^2/\text{nm}$ at 40 cm height and UV not exceed 10-4 W/m² in 180nm to 400nm.
- Effective light field >700 cm².



Touch screen panel

- Digital Hour meter showing total exposure time for current patient clearly visible by operator.
- Over temperature safety cut out is included.




- Variation in intensity over 5-6 hours < 10%.
- The irradiance ratio (min to max) is greater than 40 % on mattress.
- Green indicator light is provided to indicate that equipment is ready for normal use.
- Interruption and a restoration of the power supply does not change preset values. LED heat can be reduced by natural cooling.
- LED is protected from free fall.
- It does not topple on 10 deg inclined angle.



- Lamp life minimum 50000 hours for LED and has timer to indicate its usage.
- Each lamp unit provided Blue power LEDs emitting blue radiation between 450-480nm wavelengths for phototherapy
- Up, down and tilting of head is possible .
- The unit is mounted with castor wheels with brakes.
- The temperature of baby bed and metal surfaces does not exceed 40 deg C and 43 deg C for other accessible surfaces.
- There is intuitive method to indicate the light surface is at the appropriate treatment distance.
- Mobile stand with movable castors and height adjustment facility along with easy swivelling of source box. Unit can be used along with Infant care trolley, Radiant Warmer and Incubator.
- UP/DOWN adjustment of Over Head Unit; The phototherapy unit can be able to provide effective treatment for beds and incubators of varying heights (generally 1.0 to 1.6m). Adjustment of light intensity is provided.
- It meet IEC 60601-2-50 Medical Electrical Equipment PART 2-50: Particular Requirements for the Safety of Infant Phototherapy Equipment.

Optional Bassinet:

- Clear cabinet for observation of infant.
- Infant bassinet an integral unit which is detachable.
- Unit provide shielding of infant in the event of bulb breakage.
- Bulb mount has angle adjustment of at least 30 degrees.
- All surfaces made of corrosion resistant materials.
- Light unit tilting facility and height adjustment facility.
- The temperature of baby bed and metal surfaces does not exceeds 40 deg C and 43 deg C for other accessible surfaces.

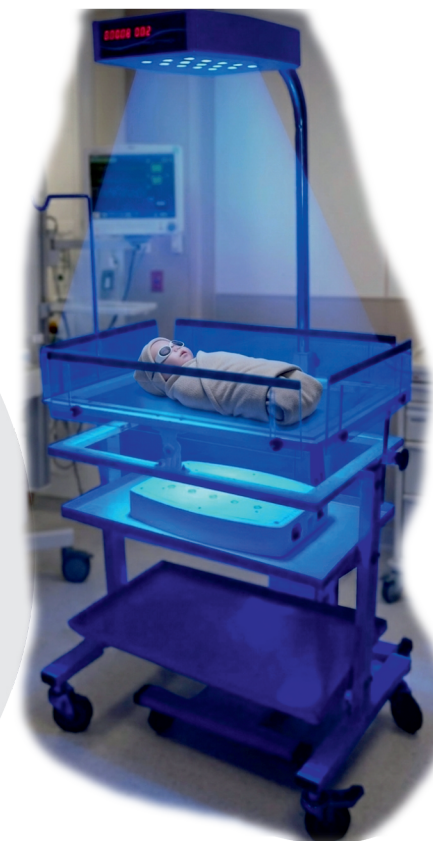



TIMER



INTENSITY
INCREASE / DECREASE BUTTON

UP/DOWN SWITCH



MEDI WAVES INC.

Works :
Plot No. 92, Sector-57, Phase-IV,
H.S.I.I.D.C. Kundli, Sonapat, Haryana (India)
Ph. : +91-8383868178

Regd. Office :
B-68, G.T. Karnal Road Industrial Area, Delhi-110033
E-mail : info@medi-waves.com

Web : www.medi-waves.com
E-mail : sales@medi-waves.com